

PATENTED

TOURNIK-IT®

MADE IN ITALY



HEMOSTASIS SYSTEM FOR MASSIVE BLEEDING



- Easy to apply even with one hand
- Adjustable to fit different limb diameters
- Effective and immediate in controlling bleeding
- Lightweight and compact, easy to carry



NSN 6515-15-036-5731
UDI DI BASE: 805759270TKY6

TOURNIK-IT It's a hemostasis system used to control severe and life-threatening bleeding. It's a simple and easy to use device especially in case of trauma to the limbs. It's works by compressing the blood vessels at the root to the limb, reducing blood flow to the affected area, thus stopping the bleeding.

THE FIRST WITH A REMOVABLE BAR

The hemostasis system can be used in extreme emergency situations, in the presence of massive bleeding; Road accidents - Gunshot or stub wounds - Serious limb injuries - Workplace accidents.

It's important to use a **TOURNIK-IT** correctly and only when necessary, as improper use can cause tissue damage or aggravate the situation.

It's essential to follow instructions and guidelines for the correct use of tourniquet and, if possible, receive specific training in the use of these devices.



Hemostatic belt



X-Pander Bar operational



X-Pander Training Bar

Permanent marker
for notes




TOURNIK-IT parts

- 1) Velcro **Hemostatic Belt** 38mm, with single buckle locking system. Stabilizing plate with "C" arch that blocks the bar, made of Poligopa, a high tenacity elastic polymer complex, resistant to UV. 2 sliding straps.
- 2) **X-Pander Bar**: removable operating bar in Polimida, a high-toughness polymer compound, resistant to UV, anatomical ipergrip design, with slot for inserting the sliding straps. Black colour.
- 3) **X-Pander Training Bar**: removable training bar in Polimida, a high-toughness polymer compound, resistant to UV, anatomical ipergrip design, with slot for inserting sliding straps. Blue colour for training.
- 4) Pocket-sized permanent **marker**.

*Innovative Italian design that uses well-known technique with the addition of an important technological solution aimed at saving and versatility. In fact when its reaches its expiration date, the removable black bar can be replaced with the blue **X-Pander Training Bar**, transforming the old no longer usable tourniquet into a training device. Identified and distinguished by the blue color of the bar, conventionally associated with training devices, thus avoiding confusion with expired black operational tourniquets.*

Additionally the X-Pander Training Bar separately to recycle damaged or the “C” arch type, from any Once the operating removed, you can easily blue bar getting distinct training can be reused

can also be purchased expired tourniquet, of manufacturer. old bar has apply the new a new and tourniquet. The bar hundreds of times.

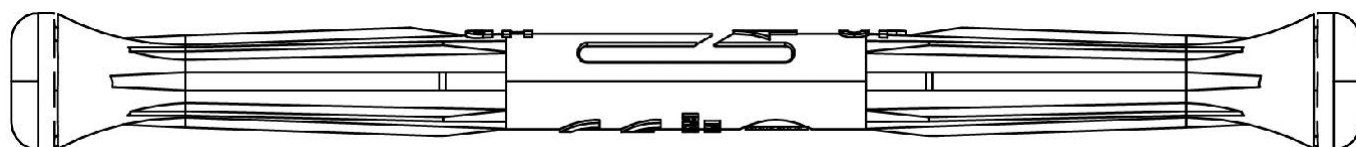

TOURNIK-IT with mounted black bar, blue bar and marker

Item code. TK-00



TOURNIK-IT with mounted black bar, marker

Item code. TK-01



TOURNIK-IT with mounted blue bar

Item code. TK-02T



Black hemostatic belt

Item code. TK-BEN



X-Pander Training Bar

Item code. XTB-01



X-Pander operating Bar



Item code. XB-01



Permanent marker

Item code. TK-PEN

TESTED IN AN INDIPENDENT LABORATORY

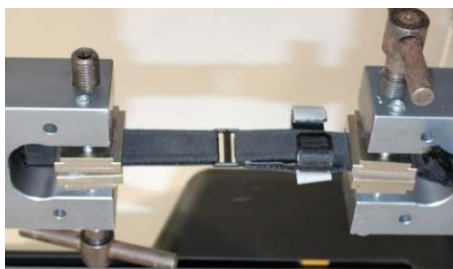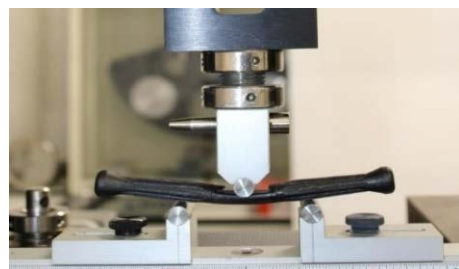

FLAMOR s.r.l.
Via Biandrate, 80
28060 SAN PIETRO MOSEZZO (NO) - ITALY
Tel. +39 0321 442441 - Fax +39 0321 392190
email: info@flamor.com
NATO NCage AE820



Compliant with regulation EU Medical Device
Regulation (MDR), 2017/745